---

# SilverCloud UF41

Non-contact Infrared Thermometer - Berührungsloses Infrarot-Thermometer - Termómetro Infrarrojo Sin Contacto - Thermomètre infrarouge sans contact - Érintés nélküli infravörös hőmérő - Termometro a infrarossi senza contatto - Bezdotykowy termometr na podczerwień - Termometru digital cu tehnologie infrarosu, non-contact

EN	User manual .....	2
DE	Benutzerhandbuch .....	5
ES	Manual de usuario .....	9
FR	Manuel utilisateur .....	14
HU	Használati utasítás .....	18
IT	Manuale utente .....	22
PL	Instrukcja obsługi .....	26
RO	Manual de utilizare .....	31

---

## General description:

- SilverCloud UF41 Non-contact Infrared Thermometer is specially designed to measure body or ambient temperature.

## Safety cautions:

- This device must only be used for the purposes described in this instruction manual.
- This device must only be used in an ambient temperature range of 5°C~40°C, the best temperature is 25°C.
- Do not expose this thermometer to extreme temperature environment of >40°C or <5°C
- Do not expose this device to electric shocks.
- Do not use this device in relative humidity higher than 80%RH.
- Do not use the device near strong electromagnetic fields such as found with cordless or cell phones.
- Keep the device away from water, heater and direct sunlight.
- Do not drop or knock the device and do not use if it is damaged.
- Do not attempt by yourself to disassemble ,repair or re design the device.

## Mode setting:

- Press “SCAN” button to activate the device and at the same time the device will start measuring the temperature. If the mode is “BODY” last time you used, the LCD maybe display “Lo”; if the mode is “SURFACE TEMP” last time you used, the LCD will show the measured result.

## Settings

- Press  button to switch between BODY and SURFACE TEMP (ambient) measuring mode.
- Press  button to switch between °C si °F
- Hold pressing  button till the screen displays “On” and you hear a confirmation beep. Release the  button. Press again  to disable the alarm.
- Press  button to switch to below functions: back light on, laser on, back light off, laser off”
- Press  button to consult memories of to consult saved data of the measured results. . The memory capacity is 32 results.

## Precaution before use:

- In order to obtain reliable and stable results, when changing the environment to a room with a significant temperature difference, it is recommended to leave the non-contact infrared thermometer to this ambient temperature for 15~20 minutes before using it.
- Depending various skin types, thickness and body position there may be different measured result. This is normal.
- Allow 1 minute interval between two measurements.

## Features:

- Alarm function once temperature is too high. The default threshold alarm value is 37.8°C (100.04°F)
- Automatically save data
- Automatically power off in 5 seconds if no operation

## Configuration:


## Display symbols:

BODY	Forehead measuring mode
SURFACE TEMP	Ambient measuring mode
°C	Centigrade
°F	Fahrenheit
	Low battery
M	Memory
	Alarm beep signal

## Taking a measurement:

- Install battery (the batteries are not included)
- For the first use or when inserting new batteries, wait 10 minutes for the warming-up of the device.
- Aim towards the forehead as the picture shows, from a distance of 5~10cm, press the scan button and measured data will be immediately displayed on the LCD. Make sure there is no air perspiration, cosmetic or cap covered on forehead. When there is sweat on forehead you can measure behind the ear lobe making sure no hair, sweater or other things cover the ear.
- Memory recall: Press “M” button to consult measured result
- Battery replacement: When LCD shows low battery symbol it means you must install new batteries. Remove the battery compartment cover and put 2AA (1.5 volt) batteries (alkaline batteries) paying attention to the polarity.

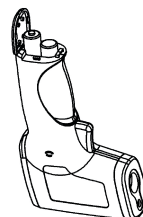---

## Trouble shooting:

If you happen to meet one of the following problems while using the device please refer to this breakdown service guide. If the problem persists please contact our customer service.

The screen display “HI” when using, which means the temperature you got is over 42°C/107.6°F on “Body” Mode.


The screen display “Lo” when using, which means the temperature you got is under 32°C/90°F on “Body” Mode.


Note: If the default measuring mode is BODY mode, when you turn on the device and measure ambient temperature you will probably see HI or Lo. In this situation, change on “SURFACE TEMP” mode. Possible causes:

1. The forehead is covered with hair or perspiration.
2. You are in an area where air currents interfere with the IR flow .
3. Temperature readings are too close together to make the body temperature have chance to boot itself.
4. The measuring distance if too far.

## Icons description:

 Measured result is over 42.9°C (Body measuring mode)

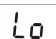 Measured result is below 32.0°C (Body measuring mode)

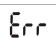 The ambient temperature is beyond the range for proper working temperature of 5.0°C - 40°C

## Specification:

Measurement range: 32.0°C - 42.9°C (BODY mode)  
0°C - 100°C (AMBIENT mode)

Accuracy ±0.3°C (32.0°C - 34.9°C)  
±0.2°C (35.0°C - 39.0°C)  
±0.3°C (39.1°C - 42.0°C)

Resolution: 0.1°C

Voice processor: English

Working temperature: 5.0°C - 40.0°C

Working humidity: Rh≤80%

Storage temperature: -25°C - 55.0°C

Storage humidity: Rh≤85%

Measure distance: 5 ~ 10cm

Power supply 3Vd.c. (2 AA alkaline batteries - not included)

Dimensions: 145 x 37 x 98 mm

Weight: 141g

## Normal temperatures according to age:

Age	Temperature °C
0-2 years	36.4-38.0
3-10 years	36.1-37.8
11-65 years	35.9-37.6
>65 years	35.8-37.5

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